
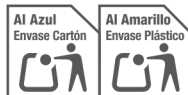
EU-A01

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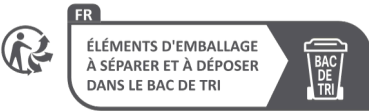1 Safety Instructions

READ & UNDERSTAND ALL INSTRUCTIONS!

WARNING: Improper use of this product may result in serious personal injury or property damage. Always follow the safety instructions in this manual for correct and safe operation.


WARNING: The use of non-approved charging cable or incompatible attachments may result in overheating, equipment damage, or personal injury.

Intended Use: This device is solely designed for oral cavity irrigation. Avoid directing water flow to sensitive areas, including the sublingual region, eyes, ears, nose, or other delicate tissues.

- This water flosser can be used by individuals with reduced physical, sensory, or mental capabilities — provided they are supervised and adequately instructed on safe usage and potential hazards.
- Place this device out of reach of infants and children. Always supervise infants and children to ensure they do not play with this product.
- Before attaching accessories, cleaning or performing maintenance, ensure the device is powered off and disconnected from the power source.
- Do not use this device in outdoor locations where aerosols (sprays) or oxygen supplies are in use.
- Remove any decorative dental accessories before operating the device.
- Avoid using iodine, saltwater, or water-insoluble concentrated essential oils with this device. These substances may impair functionality and reduce its lifespan.
- Do not attempt to repair or modify the water flosser. Unauthorized adjustments may cause injury, property damage, and void the warranty.
- Fill the tank with water below 104°F (40°C) only to prevent burns.
- If the water flosser falls into water while charging, do not touch it — unplug the power adapter immediately.
- If the power adapter is damaged, replace it before use.
- Do not press the nozzle forcefully against teeth or gums — this can cause gum injury or bleeding.
- Do not place or store the device where it can fall, slide or pulled into water.
- Do not immerse the device in water to prevent malfunction, fire, electric shock, or explosion.
- Ensure the water flosser, charging port and your hands are completely dry before charging.
- This water flosser must only be powered by the Safety Extra-Low Voltage (SELV) marked on the device.
- This product must be charged using the manufacturer-specified charging cable. Using unauthorized or incompatible chargers is prohibited.
- Before connecting to power, verify that the power adapter's voltage rating matches your local supply.
- Do not charge the device in hot or humid environments, such as bathrooms.
- Do not damage, alter, or modify the charging cable. Avoid excessive pulling, bending, or twisting. Do not place heavy objects on the cable or pinch it between objects.
- When plugging or unplugging the power cord, always hold the power adapter — not the charging cable.
- Keep the water flosser, especially the charging cable, away from heat sources.
- Place the device in a cool, dry place away from direct sunlight.
- If unused for extended periods, disconnect the charging cable to prevent fire hazards.
- Some mouthwashes contain chemicals that can weaken or crack the device's plastic parts. To prevent damage, only use clean water in the water tank.
- Before serving, disconnect the charging cable.
- For servicing, always return the product to the manufacturer, as replacement parts are not available for purchase.
- Do not operate the device if the nozzle is missing or damaged.
- Patients with periodontal disease or those who have undergone oral surgery within the last 2 months should consult a doctor before use.
- Do not replace the battery pack yourself. Replacement must be performed only by authorized technicians.
- The battery must be removed and disposed of in accordance with local e-waste guidelines. Refer to the "Disposal" section for detailed battery removal instructions.

2 Product Introduction


2.1 What's in the Box

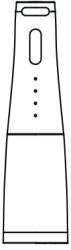Water Flosser × 1

Standard Nozzle × 1

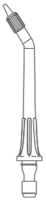Periodontal Nozzle × 1


Bubble Nozzle × 1


Tongue Cleaner Scraper × 1


Orthodontic Nozzle × 1

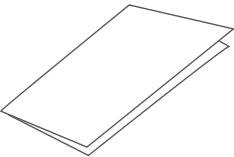User Manual x1

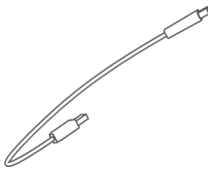
Charging Cable x1

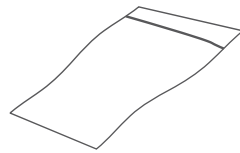Portable Bag x1

2.2 Product Overview

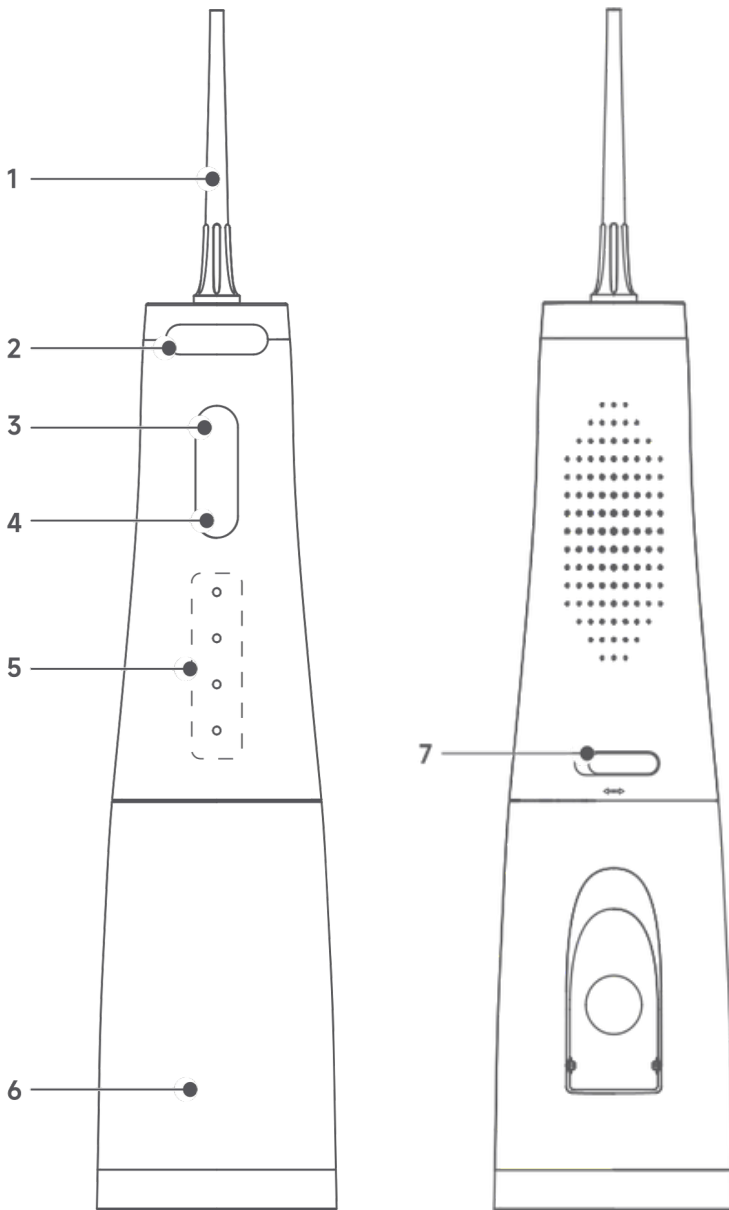- 1 Nozzle
- 2 Nozzle Release Button
- 3 Power Button
- 4 Mode Button
- 5 Mode Indicator
- 6 Water Tank
- 7 Charging Port


3 Operation

3.1 Charge the Water Flosser

Before the first use, fully charge the water flosser by following the instructions below (Fig. 1).

WARNING: Before charging, make sure the device, charging port and your hands are dry. Use only a 5V 1A power adapter (with a measured maximum output power below 15W) that complies with the GB 4706.1-2005 standard.

1. Open the charging port cover on the device. Insert the USB-A end of the charging cable firmly into the port.
2. Plug the other end of the cable into the power adapter. And insert the adapter into socket.

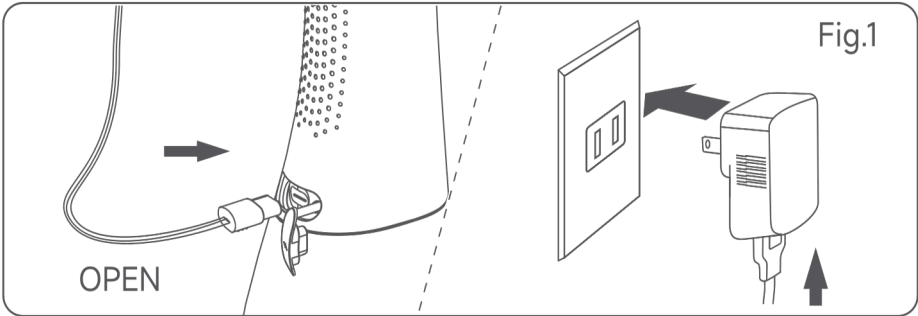NOTE:

Low Battery Indicator

- If all four indicator lights flash 3-5 times and then turn off, the device is in low battery mode and needs to be charged.

Charging Status

- While charging, the four indicator lights will flash sequentially.
- When the device is fully charged, all four indicator lights will remain solid white.


Checking Battery Level

- Press and hold the mode button to check the current battery level. The number of illuminated indicator lights shows the battery status.

3.2 Use the Water Flosser

1. Attach the Nozzle

Attach the nozzle by pushing it onto the main unit until secure (Fig. 2-1). Rotate the nozzle to adjust the direction of water flow. To replace, press the "Nozzle Release" button and repeat (Fig. 2-2).

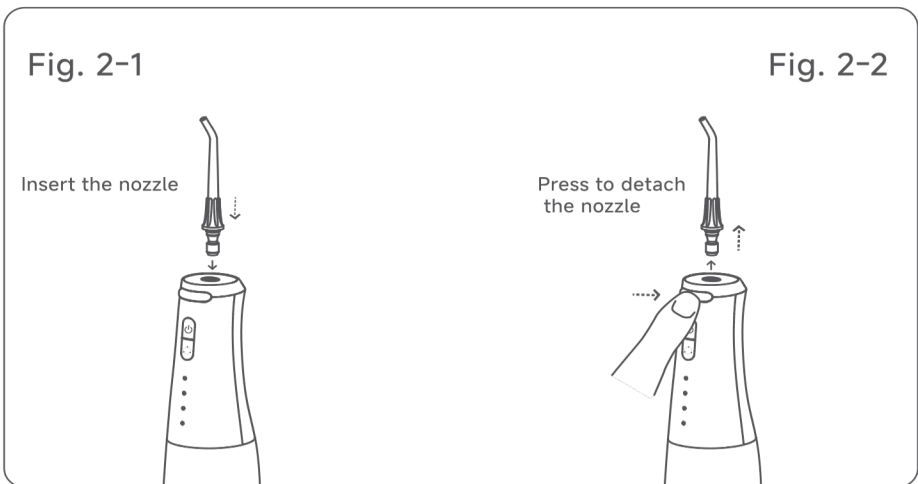2. Fill the Water Tank

a. Open the water tank cap and fill it with clean water under 40°C (104°F) (Fig. 3-1).

b. If the cap comes off, reattach it by aligning the gasket with the opening before closing securely (Fig. 3-2).

NOTE: Make sure that the device is turned off when adding water to the water tank.


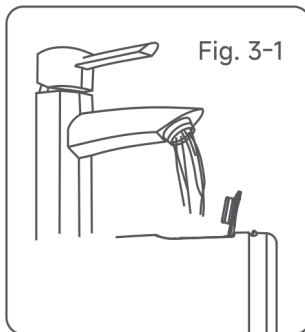
Fig. 3-2

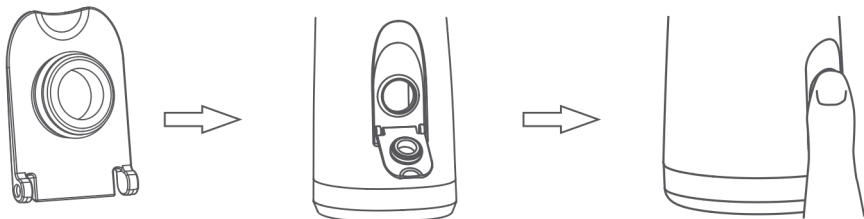3. Select Modes

Short-press the "Mode" button to cycle through flosser modes:

Low → **Medium** → **High** → **Pulse**.

Low Mode: Gentle water pressure, ideal for sensitive gums and mild cleaning.

Medium Mode: Standard water pressure, ideal for effective daily cleaning.

High Mode: Powerful water pressure, ideal for a thorough deep clean between teeth and along the gumline.



Pulse Mode: Alternates between strong and gentle pulses to deeply clean teeth and massage gums.

NOTE: For first-time users, we recommend starting with "Low Mode" for one week to adapt before switching to other modes.

NOTE: When the water flosser is turned on, it recalls the last mode used for at least one minute before the previous shutdown

4. Starting Irrigating

-When using, hold the water flosser vertically. Insert the nozzle into your mouth, align it with your teeth or gumline, and keep your mouth slightly open to allow water to flow out smoothly (Fig. 4-1, 4-2).

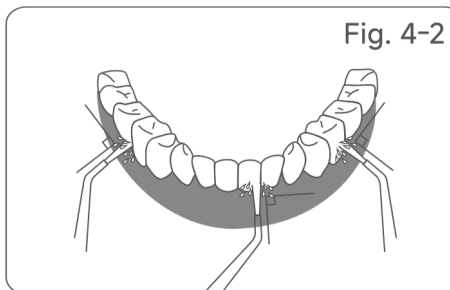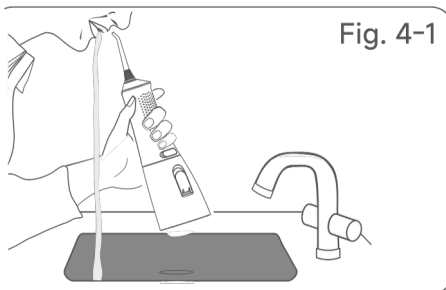- Press the "Power" button to start irrigating. Follow the instructions below during use:



a. Hold the nozzle at a 90-degree angle to your teeth and gums for optimal cleaning.

b. For best results, start with the back teeth and work systematically toward the front teeth (Fig. 5)

c. Align the nozzle with the gumline, maintaining a 90-degree angle to the gums.

d. Clean teeth, braces, crowns, and dentures thoroughly.


NOTE: Do not spray water directly into the periodontal pocket (Fig. 6).

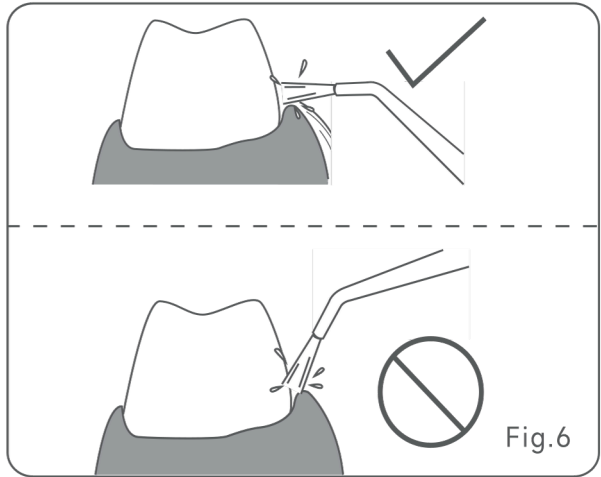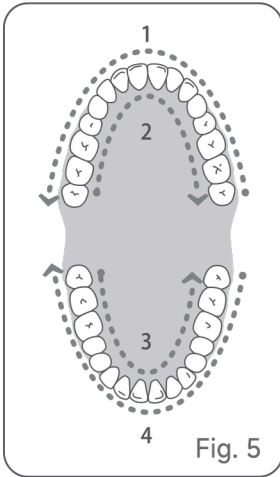- After use, press the "Power" button to turn off the device. This water flosser features a 2-minute auto-stop timer. To continue use, simply press the "Power" button again.

3.3 Post-Use Care (Fig. 7)

1. Turn off the device.
2. Pull down the water tank to detach it from the main unit.
3. Invert the tank to drain any remaining liquid.
4. Reattach the tank to the main unit.
5. Dry the Device. Wipe the exterior with a dry, lint-free cloth.

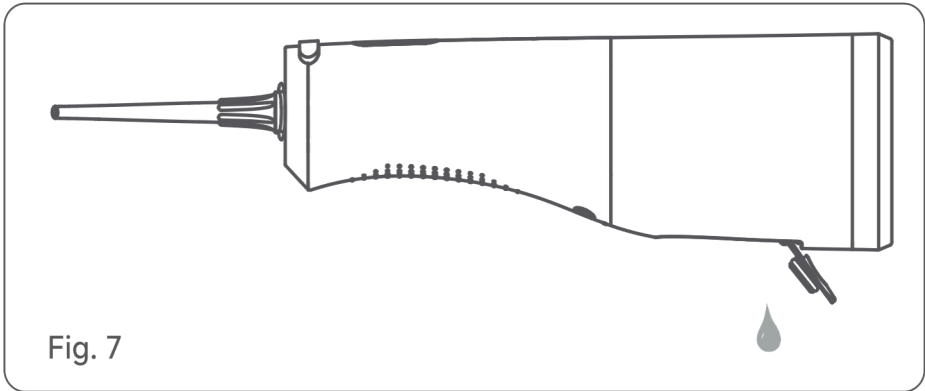NOTE: To maintain hygiene, always empty the water flosser of standing water. For long-term storage, clean and allow the device to dry completely before putting it away.

4 Maintenance

WARNING: Avoid caustic or abrasive cleaners (e.g., vinegar, lime scale removers), as they may cause damage.

WARNING: Avoid using hot water above 40°C (104°F) to clean the irrigator. Keep the device away from high temperatures and direct sunlight to ensure optimal performance.

Main Unit

1. Gently wipe away dirt using a soft cloth dampened with mild soapy water, then dry thoroughly with a clean, dry cloth.
2. Never submerge the main unit in water – moisture can damage internal components.
3. Ensure all parts are completely dry after cleaning before storage or reuse.

Nozzle

1. Rinse the nozzle thoroughly under running water and gently wipe with a soft, lint-free cloth after each use.
2. Avoid bending, pulling, or twisting the nozzle to prevent damage.
3. It is recommended to replace the nozzle every 3 months for optimal hygiene.

Water Tank

1. Rinse the tank thoroughly with clean water and gently wipe the interior and exterior with a soft, lint-free cloth. Avoid using harsh chemicals or abrasive materials.

2. Ensure the tank is completely dry before storing it for extended periods to prevent mold or bacterial growth.

Remove Hard Water

Depending on the mineral content of the water used, hard water deposits may build up on the device. Failure to remove hard water buildup will reduce product performance.

1. Fill the water tank with clean, warm water (below 40°C / 104°F).
2. Point the nozzle at the sink, turn on the device, and let it run until all the water has been drained.
3. Repeat the process by refilling the tank with clean, warm water and running the device again. Perform this cleaning every 1-3 months to maintain product performance.

NOTE: To maintain hygiene, always empty the irrigator of standing water. For long-term storage, clean and allow the device to dry completely before putting it away.

5 Disposal




For eco-friendly disposal, the device's rechargeable battery is designed for easy removal. Before discarding the product, follow the steps below to safely remove the battery. To ensure safety, always protect your eyes, hands, and fingers, and work on a stable, protected surface.

WARNING: Do not disassemble this product unless it is being permanently discarded. Unauthorized disassembly may cause irreversible damage.

WARNING: Before removing the battery, unplug the charging adapter to prevent electrical hazards.

WARNING: After removing the battery, ensure the positive and negative terminals do not short-circuit. Do not reenergize the product once the battery has been removed.

1. Remove the nozzle and water tank (Fig. 8).
2. Use a flat tool to carefully pry open the main unit bottom, and then remove the core battery component (Fig. 9).

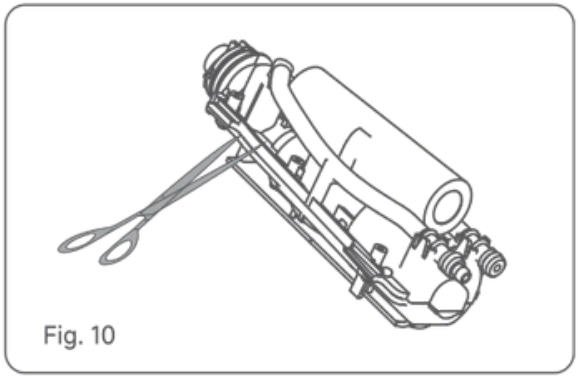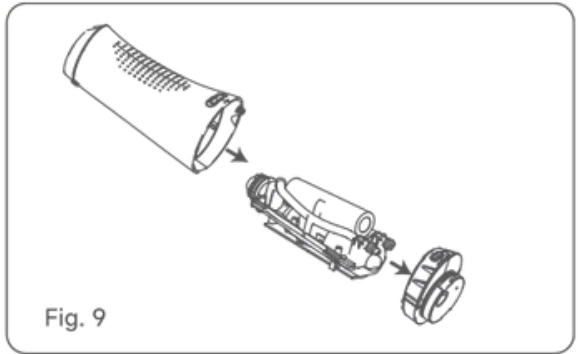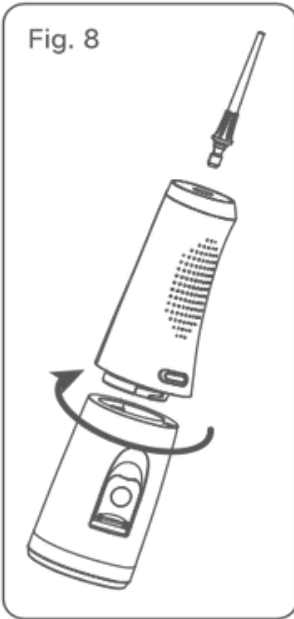3. Use scissors to cut the battery connection wires and remove the battery (Fig. 10). Recycle or dispose of the battery at a certified e-waste facility, and contact your local waste management office to confirm recycling options for the remaining device components.

6 Troubleshooting

Problem	Possible Cause	Solution
Device not working	This is your first use or the irrigator has been stored for over 3 months, a full charge is required.	Fully charge the device.
Water pressure is too low	The nozzle is deformed or is stuck.	Replace the nozzle.
	Dirt is stuck in the filter*.	Clean the filter*.

Device runs only briefly even after charging	The battery has reached the end of its service life.	Contact service center for support.
	Charging time was insufficient.	Fully charge the irrigator.
Indicator light does not turn on when charging	The charging cable may be loosely attached.	Please ensure the power is properly connected.
	The circuit is malfunctioning.	Contact service center for support.
	Low battery, insufficient to power the indicator lights.	Recharge the device for 15–20 minutes to restore normal charging indicator functionality.
No water come out	This is your first use or it has been unused for a period.	Remove the nozzle and rinse it thoroughly under running water.

*The filter is located at the lower end of the water flosser's suction hose.

7 Specification

Model	AWF14A
Product Name	Ripple 10 Water Flosser
Brand	MOVA
Rated input Voltage	5V
Rated input Power	5W 
IP-Classification	IPX7

8 Warranty

1. This product comes with 1-year free warranty starting from the date of purchase, subject to local laws and regulations.
2. During the warranty period, if the product exhibits manufacturing defects (excluding damage caused by improper use or mishandling), you may contact the MOVA after-sales service team for free repair.
3. During the warranty period, the following situations are not covered under the warranty:
 - Defects caused by normal wear and tear or normal product aging. This includes but is not limited to scratches, dents, paint peeling, or any corrosion or discoloration caused by overheating, abrasives, and chemical cleaners.
 - Products whose serial number, factory label, and other identifiers have been torn off or altered.
 - No valid purchasing document (receipt, invoice, or other) is provided to certify the date the product was purchased.
 - Damage caused by using non-specified power sources or incorrect voltage settings.
 - Damage due to accidental causes or user negligence, especially to external parts.
 - Damage resulting from improper use, storage, or maintenance by the user.
 - Damage caused by unauthorized disassembly, modification, or repair by personnel not approved by MOVA.
 - Damage caused by natural disasters (e.g., earthquakes, floods, lightning) or external incidents (e.g., fires, building collapses).

For more information and additional warranty details, please scan the QR code below to reach our support team or visit MOVA's official website: <https://www.mova.tech/>

CONTACT US
<https://global.mova-tech.com/contact>

We, Xiamen Suguang Zhiyuan Technology Co., Ltd., here by, declare that this equipment is in compliance with the applicable Directives and European Norms, and amendments. The full text of the EU declaration of conformity is available at the following internet address: www.trouver-tech.com

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmentally safe recycling. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.