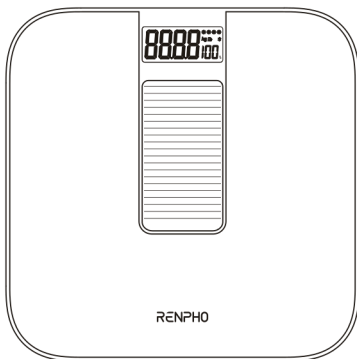
# RENPHO

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# Important Safety Instructions

## READ ALL INSTRUCTIONS CAREFULLY BEFORE USE.

1. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
2. Children should be supervised to ensure that they do not play with the appliance.
3. The scale is for household use only and not for medical or commercial purposes. It should not be used in diagnosing a disease, disorder or abnormal physical state.
4. Children under 8 or individuals who are pregnant, have osteoporosis, a pacemaker, artificial joints, portable ECG or other metal implants should use the scale for weighing purposes only and make sure to wear shoes to prevent any electric current when stepping onto the scale.
5. All body composition data provided is for reference and should not be used as a substitute for advice from a licensed healthcare professional.
6. Do not overload the scale. (Maximum Capacity: 400 lbs / 180 kg)
7. Inspect the scale for any damage before each use. Do not use a damaged scale to avoid inaccuracy or injury.
8. This device is a solar-powered scale. Charge it with direct sunlight or other light sources. It does not support USB charging.
9. Turn on the power switch before the first use. Switch the scale off if it is not used for a long time.
10. To activate the scale, step on it without shoes.
11. Do not use the scale in the environment with sudden temperature changes.
12. To ensure accuracy, use the scale on a hard, flat and dry surface.
13. To avoid interference or danger, do not use the scale in an environment near equipment with strong electromagnetic fields.
14. To prevent injury, do not stand on the edges of the scale, jump on or stomp it.

15. The scale is a delicate and precise measuring device. Handle with care to avoid breakage.
16. Place/use the scale in a well-lit area (At least 2 " ☀ " appear on the display.) to ensure it receives sufficient light for prolonged use.
17. Place the scale in the area where 4 " ☀ " appear on the display to obtain faster charging when it is in low battery.
18. Avoid splashing water on the solar panel on the scale.
19. The scale is not waterproof. Keep it dry and never submerge it in water. To clean the surface, use a damp cloth or glass cleaner. Do not use soap or chemical cleaning agents.
20. Do not attempt to disassemble the scale for any reason, as this may cause damage or void the warranty.

### **About the Battery**

1. Do not expose the appliance or battery to excessive temperatures.
2. Be aware of the risk of terminals of the battery-operated appliance or battery being short- circuited by metal objects.
3. This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of.
4. Do not dismantle, open or shred battery.
5. Do not subject cells or batteries to mechanical shock.
6. In the event of a cell leaking, do not allow the liquid to come in contact with the skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical advice.
7. Secondary cells and batteries need to be charged before use. Always use the correct charger and refer to the manufacturer's instructions or product manual for proper charging instructions.
8. Do not leave a battery on prolonged charge when not in use.
9. After extended periods of storage, it may be necessary to charge and discharge the cells or batteries several times to obtain maximum performance.
10. CAUTION - The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not recharge, disassemble, heat above 100°C (212°F).

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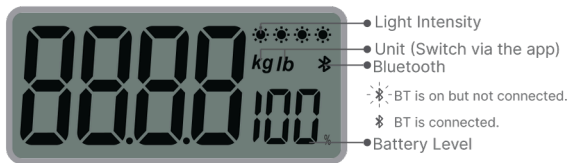
Product Weight: 1.3kg / 2.9lb

## 2. Display

- To ensure the scale receives sufficient light for prolonged use, place/use it in a well-lit area. (At least 2 "☀️" appear on the display.)
- To obtain faster power storage, place the scale in the area with direct sunlight (4 "☀️" appear on the display.) when it is in low battery.

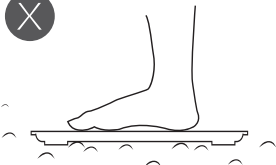
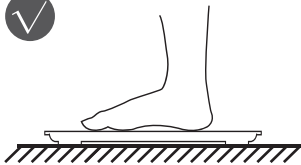
### 3. Icons

 Overloaded (Max. 180kg/400lb)	 Low battery (Place the scale in the sun timely.)	 Measuring body fat	 Calibrating
---	---	---	--

## Using the Product

### 1. Before Use

- a. Place the scale on a dry, flat, solid and even surface. Avoid any wet, slippery, carpeting or uneven flooring.

- b. Turn ON the power switch on the back of the scale.  
c. Calibrate: Step on the scale platform gently and then step off. The scale's display will show "CAL" and then "0.00", indicating that calibrating is complete.  
d. Start measuring.

## 2. Measuring

- Step on the scale with dry and bare feet and wait for 3-5 seconds.
  - When the number shown on the scale's display flashes twice, measuring is complete.
- \* The measurements can be viewed via the Renpho Health App if the connection is successful.

### Tips:

- To ensure accuracy, please try to measure at the same time and under the same conditions.
- The scale will be on the standby status if left idle for 10 seconds. Step on the scale gently to activate it.

## App Guides

### 1. Download the Renpho Health App

Scan the QR code below or search for "Renpho Health" from the App Store or Google Play to download the app.

System Requirements: iOS 12.0 or later, Android 7.0 or later.

Renpho Health

Renpho Health

Note: Due to continuous updates and improvements, the "Renpho Health" App may appear slightly different.

## 2. Log in/Register

- Log in directly with an existing Renpho App account.
- If you don't have any Renpho App account, tap **"Sign up"**. Register a new account using your email address and complete your profile information, which is required for calculating the measurements only.

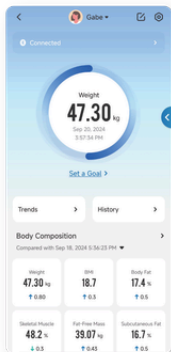
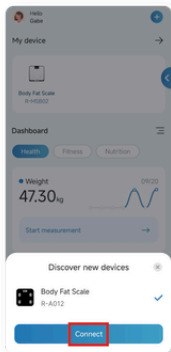
**Note:** You can also log in using Guest Mode if you do not wish to register an account. However, your measurements will not be saved if you log in on another device.

## 3. Add Device

- \* Enable Bluetooth on your smartphone before adding a device. (Location is required by Android System.)




### Method 1:

- Step on the scale and a window will pop up.
- Tap **"Connect"** to pair your smartphone with the scale.



## Method 2:

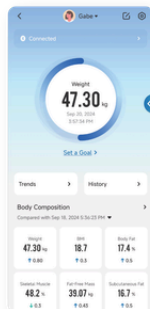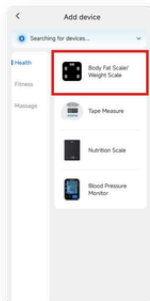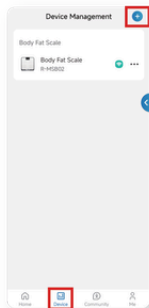- Tap “+” on the top right corner of the “Device Page” to connect manually.
- Tap “Body Fat Scale”.
- Step on the scale following the in-app instructions until a window pops up on the screen.
- Tap “+” to connect. (The Bluetooth icon will turn solid on the scale’s display once the pairing is successful.)

## 4. Switch Units

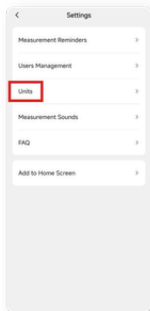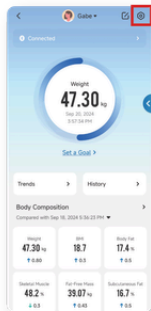To switch units, tap “⚙️” and then tap “Units” to choose the unit according to your needs.

## FAQ

### 1. How can I turn on the scale?

Turn ON the power switch on the back of the scale and step onto activate it without shoes.

### 2. Do I need to replace the battery?

No, you don't. The scale is powered by solar power and other light sources.

### 3. How can I charge the scale?

Please place the scale where light is noticeable. The stronger the light intensity, the faster the charge.

### 4. How can I know the current battery level?

Check the battery percentage on the scale's display. 100% means that the scale is at full charge.

### 5. How can I know the current light intensity?

Check the number of "  " on the scale's display. The more "  " there are, the stronger the light intensity is.

### 6. Why does the scale stop working?

The scale may have run out of power. Please place the scale in a well-lit place for charging.

## Warranty Policy

Your RENPHO product purchase is covered by a one year limited manufacturer warranty from the date of delivery. For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty. If you choose not to register your product, it will not void the product warranty.

## Customer Service

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

 TEL: +44(0)7434 666088 (UK ONLY)  
Monday-Friday 9:00AM-5:00PM (GMT)  
+1(844) 417 0149 (US&CA)  
Monday-Friday 9:00AM-4:30PM  
1800 MY RENPHO (1800 69 7367) (AU ONLY)  
Monday-Friday 9:00AM-6:00PM (AEST)  
+49 40 23969285 (DE ONLY)  
Monday-Friday 9:00AM-6:00PM (CET)

 Email: [support-eu@renpho.com](mailto:support-eu@renpho.com) (UK&EU)  
[support@renpho.com](mailto:support@renpho.com) (US&CA)  
[support-au@renpho.com](mailto:support-au@renpho.com) (AU)

\*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

### **FCC Regulatory Compliance**

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

### **ISED Regulatory Compliance**

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

1. This device may not cause interference.
2. This device must accept any interference, including interference that may cause undesired operation of the device.

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1. Meaning of the crossed-out wheeled bin symbol  
The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.
2. Options for returning old appliances  
Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.
3. Personal data  
All end-users of WEEE are reminded that you are responsible for deleting personal data from the WEEE to be disposed of.

1. The symbol of the crossed-out garbage can on batteries or accumulators indicates that they must not be disposed of in household waste at the end of their life.
2. You are legally obliged to return old batteries and accumulators after use. You can do this free of charge at a collection point in your area. Addresses of suitable collection points can be obtained from your city or local government.
3. Batteries may contain substances that are harmful to the environment and human health. The separate collection and recycling of old batteries and accumulators is intended to avoid negative effects on the environment and human health.
4. Please avoid the generation of waste from old batteries as far as possible, e.g. Please avoid littering public spaces by not carelessly leaving batteries or electrical and electronic equipment containing batteries.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/53/EU, Directive 2011/65/EU, Radio Equipment Regulations 2017, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address: <https://renpho.eu/pages/compliance-data>

Maximum RF power: -3.06 dBm  
Operation frequency: 2402 MHz to 2480 MHz

# RENPHO