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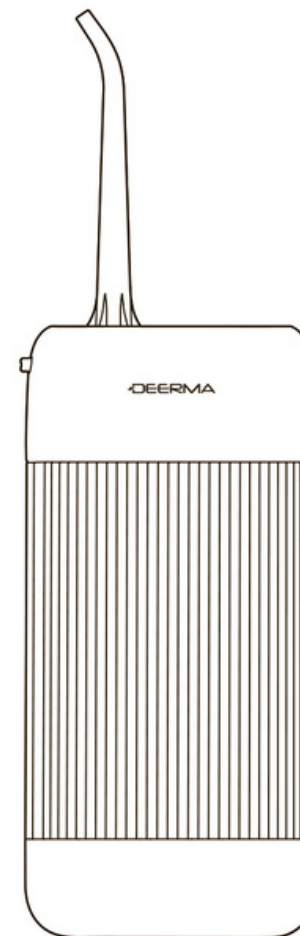
### Safety Instructions

1. If you are currently undergoing oral treatment or suffer from severe periodontal disease, consult your dentist before use.
2. If gum bleeding persists, it may indicate a gum condition such as periodontal disease; stop using the product and consult a dentist.
3. This product is a personal-care item and does not replace medical devices; do not share it with others.
4. Never operate the device without the nozzle attached.
5. Do not use water hotter than 40 °C or colder than 10 °C in the tank for oral cleaning to prevent scalding or frostbite.
6. Do not charge in humid environments.
7. Do not allow foreign objects to fall into or be inserted into the charging port, as this may cause electric shock or a short circuit.
8. Ensure both the charging port and the charging cable are dry before charging; otherwise, a short-circuit hazard may occur.
9. Do not use this product for purposes other than oral cleaning.
10. Do not press the nozzle forcefully against gums or teeth to avoid gum injury.
11. Never disassemble, modify, or forcibly bend the product yourself.
12. This product contains a built-in rechargeable battery; do not throw it into a fire, expose it to heat, or charge, use, or store it at high temperatures, as this may cause overheating, fire, or explosion.
13. This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
14. Children should be supervised to ensure that they do not play with the appliance. Children being supervised not to play with the appliance.
15. This appliance contains batteries that are only replaceable by skilled persons.
16. Must be only supplied at safety extra low voltage corresponding to the marking on the appliance.

### Care & Maintenance

#### 1. Daily use

- After each use, always empty the remaining water from the tank, then run the unit dry for a few seconds to expel any residual moisture and prevent bacteria, odors, or malfunctions.
- Rinse the nozzle under running water to remove food debris and deposits.

## 2. Long-term non-use

- Completely empty and dry the tank, water lines, and nozzle.
- Remove the nozzle and store it separately to prevent deformation.

## 3. Battery care

- Ensure the battery charge is between 40 % and 60 % before storage to avoid deep-discharge damage.
- Recharge at least once every 3 months.

## FAQ

### 1. Cannot power on

Check that the product has sufficient battery level.

### 2. No water spray after powering on

The nozzle may be clogged—clean or replace it.  
Ensure the water tank contains enough water.

### 3. Water tank leaks

Ensure the water tank is correctly installed.

All products bearing this symbol are waste electrical and electronic equipment (WEEE as in directive 2012/19/EU) which should not be mixed with unsorted household waste. Instead, you should protect human health and the environment by handing over your waste equipment to a designated collection point for the recycling of waste electrical and electronic equipment, appointed by government or local authorities. Correct disposal and recycling will help prevent potential negative consequences to the environment and human health. Please contact the installer or local authorities for more information about the location as well as terms and conditions of such collection points

## Product parameters

Product name: Portable water flosser

Model: CY01

Operating voltage: 3.7 V

Operating power: 8 W

Charging voltage: 5 V =

Charging current: ≤0.5 A

Battery capacity: 1500 mAh

Product size: 140×66×52 mm