ARTEMIS

BLUSH EDITION

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Thank you for your purchase.
Please read this manual carefully before use and follow all instructions. As we continually improve our products, specifications may change – always refer to your actual device if features differ from this manual. Retain this guide for future reference.

IMPORTANT SAFETY INSTRUCTIONS

READ ALL INSTRUCTIONS CAREFULLY BEFORE USING.

1. The provided eye cups enhance safety by blocking infrared light and other light sources. If using the device without the eye cups installed, avoid looking directly at the light source, limit prolonged exposure, and keep your eyes closed during use.
2. When using for the first time, please fully charge the controller.
3. The mask contains circuits and LED modules—avoid squeezing, folding, or placing heavy objects on it.
4. Use only a power supply with $\leq 15\text{W}$ output.
5. Discontinue use if discomfort occurs.
6. For better absorption, use with disposable facial masks containing natural ingredients or skincare products like acne/whitening treatments to help enhance the effects.
7. The mask contains silicone and may have a slight odor when first opened—this is normal and will disappear quickly.

WARNINGS

- Do not leave the device unattended when turned on to avoid the risk of fire or burns.
- This device is not designed for use by persons with limited physical, sensory, or mental abilities, or lack of experience and knowledge, unless they are trained, supervised, or used together to avoid fire hazards or burns.
- This device is not intended for use by children. Children should be supervised to ensure that they do not play with the device to avoid the risk of fire and burns.
- Do NOT wash electrical parts with water or other liquids.
- Water and electricity are a dangerous combination, and in order to avoid the danger of the exit point:
 - Do not use the device in a humid environment (e.g., in a bathroom or shower or pool attachment).
 - Do not allow water to enter the device.
- Do not use the device if it appears to be operating incorrectly, and contact customer support.
- Do not use the device in combination with any medications which may result in damage to skin tissue.
- Do not expose any parts of the device to chemical solvents or harsh cleaning fluids.

- The device contains silicon and plastics that come into contact with the skin during use. Do not use the device if you have a known allergic reaction to contact with silicone or plastics.
- Avoid fluids on the light panel and controller. Discontinue use if there is water in the device.
- Opening or altering the device may damage the operating circuits and will result in voiding the warranty.
- Please keep the unit out of reach of children, since the strap can cause strangulation.
- Use of controls or adjustments or performance of procedures other than those specified herein may result in HAZARDOUS radiation exposure.

This device is for cosmetic purposes only and is not intended to diagnose, treat, cure, or prevent any disease. It is not a medical-grade LED therapy device and has no therapeutic or clinical effects.

About the Battery

1. Do not expose the device or battery to excessive temperatures.
2. Be aware of the risk of terminals of the battery-operated device or battery being short-circuited by metal objects.
3. This device contains batteries that are non-replaceable. When the battery is at end of life, the device shall be properly disposed of.
4. Do not dismantle, open or shred battery.
5. Do not subject cells or batteries to mechanical shock.
6. In the event of a cell leaking, do not allow the liquid to come in contact with the skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical advice.
7. Secondary cells and batteries need to be charged before use. Always use the correct charger and refer to the manufacturer's instructions or product manual for proper charging instructions.
8. Do not leave a battery on prolonged charge when not in use.
9. After extended periods of storage, it may be necessary to charge and discharge the cells or batteries several times to obtain maximum performance.
10. CAUTION - The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not recharge, disassemble, heat above 100°C (212°F).

Battery Producer Labelling

Manufacturer: Dongguan Leijie Energy Technology Co., Ltd.

Address: Room 205, No. 8 Shangheng Road, Chang'an Town, Dongguan City, Guangdong Province


Mail: 2933498368@qq.com

Website: www.hdbattery.cn

SAVE THESE INSTRUCTIONS

CONTAINS

① LED silicone mask

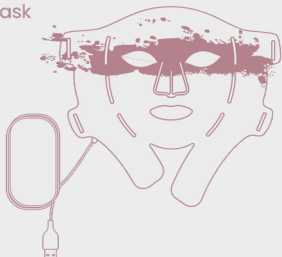
② Eye cup

③ USB charging cable

④ Controller

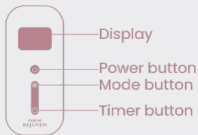⑤ Adjustable head strap

PREPARATION

Fully charge the device before first use (about 2.5 hrs) - indicators flash during charging and stay solid when fully charged (Note: Cannot be used while charging).

Please follow the steps in order

Step 1

Clean and dry your face.

Step 2

Properly install eye cups and head straps.

Step 3

Connect the controller.

Step 4

Wear the mask, adjust the straps for a secure fit.

OPERATION

Light Sensitivity Skin Test

Before first use, perform the following test:

Step 1

Place the device on your arm and turn it on for 10 minutes. If redness or itchiness occurs, you may have photosensitive skin. Discontinue use or consult a dermatologist.

Step 2

Clean and dry your face (remove all makeup) to avoid interference from light-sensitive cosmetic ingredients. Apply the mask to your face.

Step 3

Check for any adverse reaction after use, mild redness is normal and should fade within 1–2 hours. If burning or discomfort occurs, stop immediately and seek medical advice.

1. Power On/Off

Turn On: Press and hold the Power button.

Default mode: Red + NIR Light (Intensity Level 1) with a 10-minute timer.

Turn Off: Press and hold the Power button again

2. Adjust Intensity

Press the Power button to cycle through 5 intensity levels.

3. Switch Light Modes

Press the Mode button to toggle between Red (R), Blue (B), Green (G), and Infrared Light (NIR) combinations:

R + NIR Mode (630nm + 850nm)

B + NIR Mode (415nm + 850nm)

G + NIR Mode (520nm + 850nm)

R + B + G + NIR Mode (630nm+ 415nm+ 520nm + 850nm)

4. Set Timer

Press the Time button to switch between 10, 20, or 30-minute sessions.

Frequency

Initial Phase: 3-4 sessions/week (10 minutes each).

Maintenance Phase: 1-2 sessions/week after visible results.

SPECIFICATION

Product Name: LED Light Soft Mask

Model: E100B-1

Input: 5V $\overline{=}$ 1A

Battery Voltage: DC 3.7V

Battery Capacity: 2000mAh 7.4Wh

Timer: 10/20/30 min (10 min recommended)

Charging Time: 2.5H

Operating Time: >45 min

Recommended Ambient Temperature: 50°F - 104°F (10°C - 40°C)

CLEANING

1. Do not soak the mask in water. Gently wipe the surface with a damp cloth.
2. Remove head straps for separate washing.
3. Never use corrosive solvents (alcohol, bleach, etc.) to clean the mask.

STORAGE

1. After your session, store the mask in its box (or a cool, dry place).
2. Avoid direct sunlight or high temperatures to prevent material degradation.

TROUBLESHOOTING

Q: Why does the device shut off during use?

- Automatic shutdown occurs when session time ends.
- Check battery level.

Q: Why is the light blinking during use?

- Ensure the mask is securely connected to the controller.

Q: Which intensity level should I choose?

- Start at Level 1 (lowest) and gradually increase based on skin comfort.

Q: Can I use it daily for faster results?

- No. Follow the recommended schedule. Overuse may irritate skin without improving efficacy.

Q: Why am I not seeing the same results as others?

- Results vary due to individual skin thickness and consistency of use.

WARRANTY (1-YEAR-PERIOD)

The manufacturer's warranty is void under the following circumstances but not limited to:

- Improper or inadequate maintenance or modification.
- Accident, misuse, abuse, contamination, or other external causes.
- The use of unspecified adapter and accessories.
- Loss or damage in transit.
- Damages that occur as a result of one's failure to follow the instructions.

- This warranty does not apply to expendable or consumable parts and does not extend to any product from which the serial number has been removed.

CUSTOMER SERVICE

We have an award-winning warranty, exchange, and customer service team that guarantees hassle-free solutions to any issue you might have within business hours.

TEL: +1(844) 417 0149 (US&CA)

Monday-Friday 9:00AM-4:30PM

1800 69 7367 (AU ONLY)

Monday-Friday 9:00AM-6:00PM (AEST)

+44(0)7434 666088 (UK ONLY)

Monday-Friday 9:00AM-5:00PM (GMT)

E-MAIL: support@renpho.com

* For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

1. The symbol of the crossed-out garbage can on batteries or accumulators indicates that they must not be disposed of in household waste at the end of their life.
2. You are legally obliged to return old batteries and accumulators after use. You can do this free of charge at a collection point in your area. Addresses of suitable collection points can be obtained from your city or local government.
3. Batteries may contain substances that are harmful to the environment and human health. The separate collection and recycling of old batteries and accumulators is intended to avoid negative effects on the environment and human health.
4. Please avoid the generation of waste from old batteries as far as possible, e.g. Please avoid littering public spaces by not carelessly leaving batteries or electrical and electronic equipment containing batteries.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/30/EU, Directive 2011/65/EU, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address: <https://renpho.eu/pages/compliance-data>

RENPHO
REJUVEN