
RENPHO

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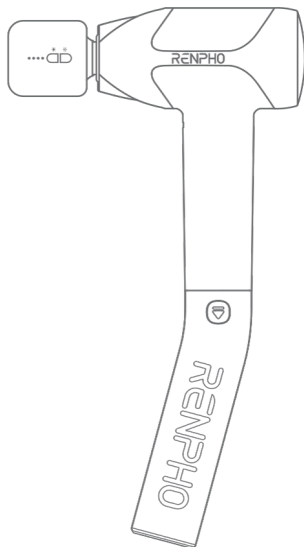Thank you for choosing the RENPHO Massage Gun! Ergonomically designed with a detachable extension handle for hard-to-reach areas.

Please carefully read this user manual and use the device only as directed. RENPHO assumes no responsibility for any injuries caused by misuse or negligence in following the instructions provided.

Please keep this manual for future reference.

Important Safety Instructions

READ ALL INSTRUCTIONS CAREFULLY BEFORE USING.

When using this device, the following basic precautions should be always followed:

1. The device is intended for household use only and not intended for medical use. It should not be used in treatment or alleviation of disease or abnormal physical state.
2. This device is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the device by a person responsible for their safety.
3. Children should be supervised to ensure that they do not play with the device.

DANGER - To reduce the risk of electric shock:

1. Do not attempt to retrieve the device if it falls into water.
2. Do not use while bathing or in a shower.
3. Do not use the device with wet hands.
4. Do not place or store the device in areas where it may fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.

WARNING - To reduce the risk of fire, electric shock, or injury to persons:

1. Please consult your doctor before using if you are a user of

the following persons:

- 1) Children under 8
- 2) Pregnant women
- 3) Persons with blood stasis
- 4) Persons with fractures
- 5) Persons with heart disease, a pacemaker, defibrillator, or metal contents in the body
- 6) If you have experienced joint surgery within 90 days, do not use the device on your body within 3 inches of the surgical site.

If you feel unwell or have any unusual discomfort, please stop using the device immediately.

2. Do not use the massage gun on

- 1) Sensitive skin, swollen or inflamed areas, in areas where there is poor circulation, where skin eruptions are present, or in the area of unexplained abdominal pain
 - 2) Head, face, cervical vertebrae, body bones, varicose veins, tumor site, throat area, any possible fracture site, positions close to the bony areas: the humerus, the back of the foot, the back of the hand, etc
 - 3) Genital area (Do not use this device for sexual purposes. It is a powerful massaging device and may cause severe bodily harm.)
 - 4) Unconscious or sleeping persons
-
3. Do not operate under a blanket or pillow. Excessive heating may occur and cause fire, electric shock, or injury to persons.
 4. Never operate the device with the air openings blocked. Keep the air openings free of lint, hair, etc.
 5. Close supervision is necessary when this device is used by, on or near children, invalids, and disabled persons.
 6. Use this device only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
 7. Never operate this device if it is not working properly, has been damaged, or dropped into water.
 8. Never drop or insert any object into any opening.
 9. Do not operate in the areas where aerosol (spray) products are being used or where oxygen is being administered.

10. Do not keep the device near a magnetic field or in the areas where it can be affected by a magnetic field as it may cause batteries to discharge.
11. Avoid touching the gap between the massage head and body with your finger or hair to avoid pinching.
12. Keep the device away from heat and avoid using it in places where the surrounding temperature is relatively hot.
13. Do not drop, apply excessive force or heavy weight over the device to avoid malfunction and/or damage.
14. If a high temperature or foul odor occurs while using the device, stop using the device immediately.
15. Do not use the device in close proximity to long hair, loose clothing or jewelry.
16. Avoid applying excessive pressure on the massage head and refrain from massaging the same area continuously for an extended period of time.
17. The device has a heated surface. Persons insensitive to heat must be careful when using the device.

About the Battery

1. Do not expose the device or battery to excessive temperatures.
2. Be aware of the risk of terminals of the battery-operated device or battery being short-circuited by metal objects.
3. This device contains batteries that are non-replaceable. When the battery is at end of life, the device shall be properly disposed of.
4. Do not dismantle, open or shred battery.
5. Do not subject cells or batteries to mechanical shock.
6. In the event of a cell leaking, do not allow the liquid to come in contact with the skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical advice.
- 7 Secondary cells and batteries need to be charged before use. Always use the correct charger and refer to the manufacturer's instructions or product manual for proper charging instructions.
8. Do not leave a battery on prolonged charge when not in use.
9. After extended periods of storage, it may be necessary to charge and discharge the cells or batteries several times to obtain maximum performance.

10. CAUTION - The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not recharge, disassemble, heat above 100°C (212°F).

SAVE THESE INSTRUCTIONS

Battery Producer Labelling

Manufacturer: SHEN ZHEN JIAJINYUAN TECHNOLOGY CO., LTD

Address: Building A3, Fuying Second Industrial Zone, Jian'an Road, Zhancheng, Fuhai, Bao'an, Shenzhen, Guangdong, P.R. China

Website: www.zgszjy.com

Email: wsh@gdszjy.com

About the Product

1. What's in the Box

1 x Massage Gun

6 x Massage Head

1 x Extension Handle

1 x Travel Case

1 x USB-C Charging Cable

2 x Rubber Sleeve

2. Product Parameter

Model: R-WMG02

Input: 5V --- 2A

Battery Voltage: DC 7.4V

Speed: 1600-2800r/min

Battery Capacity: 2500mAh 18.5Wh

Auto-off Timer: 10 minutes

Charging Time: about 2.5 hrs

Operating Time: about 3 hrs

Noise: $\leq 55\text{dB}$

Dimension: 14.3 x 5.5 x 2.7in /

364 x 140 x 68mm (with extension handle)

7.7 x 5.5 x 2.7in /

195 x 140 x 68mm (without extension handle)

Weight: 1.3lb / 608g (massage gun only)

Heating & Cooling Head

Input: 5V --- 2A

Battery Capacity: 700mAh 2.59Wh


Auto-off Timer: 10 min (heating); 5 min (cooling)

Charging Port: USB-C

Charging Time: about 2 hrs

Operating Time: about 1.5 hrs


3. Product Overview

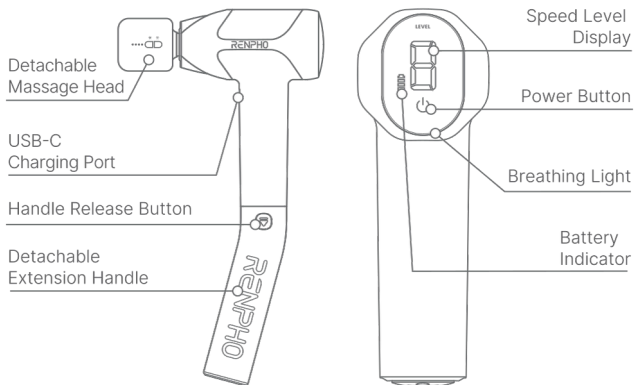Before Using

1. Charge the Device

* To ensure optimal performance, please fully charge the device before its first use.

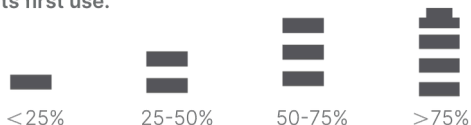Indication of remaining battery power

- Please charge the device using the provided USB-C charging cable and a DC 5V 2A charger (not included) when the last battery indicator flashes.
- It takes about 2.5 hours to fully charge the device. When charging the device, the battery indicators will display as below.


<25% ■■■■■	4 indicators flash
25%-50% ■■■■■	1 indicator on, 3 indicators flash
50%-75% ■■■■■	2 indicators on, 2 indicators flash
75%-85% ■■■■■	3 indicators on, 1 indicator flashes
>85% ■■■■■	4 indicators on

Note:


- For safety reasons, the device cannot be turned on while charging.
- The device runs on a single charge for 3 hours (depending on speed level and pressure applied during use).

2. Attach/Detach the Extension Handle

- a. To attach the extension handle, push it in until you hear a "click", indicating it is installed correctly and securely. Ensure it is installed in the correct direction as shown below.

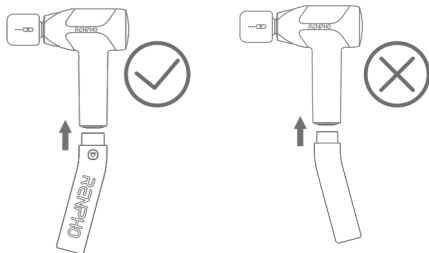- b. To detach the extension handle, press and hold the handle release button and pull the extension handle out, as shown below.

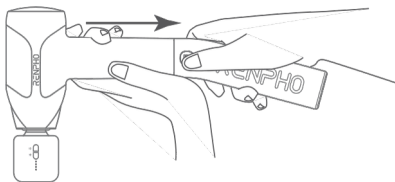3. Choose an appropriate massage head and attach it.

Ensure the massage head is fully inserted into the slot.

Note: Do not install or replace the massage head while the device is operating.

Heating & Cooling massage head

(Rechargeable)
For all body parts.


Air cushioning head

For sensitive muscles.


Ball head

For large muscle groups like quads, glutes, extreme muscles, and joint soreness.

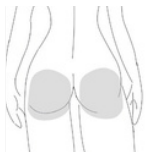U-shaped head

For the spine, neck and Achilles tendon.

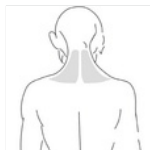Flat head

For the shoulder blades and IT bands or for a relaxing full-body massage.

Bullet head

For joints, deep tissues, trigger points and small muscle areas like feet and wrists.

Rubber sleeve

Each massage head is pre-installed with a rubber sleeve that ensures the massage head is securely attached to the massage gun. If it is damaged, replace it with a spare one.

Using the Product

1. **Press and hold** the Power button for 2s to turn on the device, speed level 1 will be activated by default.


2. **Press** the Power button to cycle through the speed levels among Level 1/2/3/4/5/Au-to-speed mode.

Level 1: 1600r/min Level 2: 1900r/min

Level 3: 2200r/min Level 4: 2500r/min

Level 5: 2800r/min

Auto speed mode: 1600-2800r/min (Automatically speed up and down based on the amount of pressure applied.)

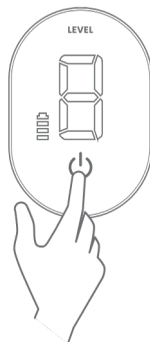Breathing Light

Status	Pressure Level	Indication
Blue	Gentle	Normal operation
Orange	Medium	Stall warning
Red	Strong	Stalled

Working Principle: The higher the speed level, the higher the stall force (Stall force is the amount of pressure the device can withstand without stalling or stopping.)

3. To turn off the device, **press and hold** the Power button for 2s.

Auto-off Timer: The device will automatically shut off after 10 minutes of operation. It is not recommended to use the device for more than 30 minutes at a time.

Tips:

- a. Massage stiff muscles for up to 2 minutes to promote relaxation.
- b. It is recommended to start with the lowest setting and gradually increase the speed as you adapt to the intensity.
- c. Do not apply too much pressure on the device during your session, or the device will stall to prevent injury or bruising.

Heating & Cooling Massage Head

- a. Press and hold to enable the heating function. RED indicator will turn on.
- b. Press to cycle through heating levels: 1-2-3-4.
- c. Press and hold again to disable the heating function.
Heating Temperature: 98.6-113°F/37-45°C. Timer: 10 min

- a. Press and hold to enable the cooling function. BLUE indicator will turn on.
- b. Press to cycle through cooling levels: 1-2-3-4.
- c. Press and hold again to disable the cooling function.
Cooling Temperature: 46.4-60.8°F/8-16°C. Timer: 5 min

Cleaning and Maintenance

1. Place the device in a cool, dry place when not in use.
2. Always turn off the device before cleaning. Wipe with a soft cloth slightly dampened with non-abrasive cleaners. Keep away from all solvents and harsh detergents.
3. Do not attempt to disassemble or repair the device. The device does not contain user-serviceable parts.

Troubleshooting Guide

Problem	Cause	Solution
The device cannot be turned on.	Press the Power button.	Press and hold the Power button for 2s to turn it on.
	The battery runs out.	Charge the device.
The device shuts down suddenly.	Auto-off timer: The device will automatically shut off after 10 minutes of operation.	Restart the device for another session.
	Built-in locked-motor protection: When the massage head is pressed with excessive force, the protection will be triggered, shutting the device off.	
	Built-in overheat protection: It will automatically shut off when overheating occurs.	Wait for the device to cool down before turning it on again.
The device fails to work.	The massage head is inserted too tightly.	Pull out the massage head a little and turn on the device.
Fail to attach the extension handle.	The extension handle is not oriented correctly.	Check the instruction on page 7.
Fail to detach the extension handle.	The handle release button is not pressed.	Press and hold the button and pull the extension handle out.
If you encounter any other problems, please feel free to contact our customer service.		

Warranty Policy

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty.

If you choose not to register your product, it will not void the product warranty.

Customer Service

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

 **TEL: +44(0)7434 666088 (UK ONLY)**
Monday-Friday 9:00AM-5:00PM (GMT)
+1(844) 417 0149 (US&CA)
Monday-Friday 9:00AM-4:30PM
1800 MY RENPHO (1800 69 7367) (AU ONLY)
Monday-Friday 9:00AM-6:00PM (AEST)

 **Email: support-eu@renpho.com (UK&EU)**
support@renpho.com (US&CA)
support-au@renpho.com (AU)

*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

1. The symbol of the crossed-out garbage can on batteries or accumulators indicates that they must not be disposed of in household waste at the end of their life.
2. You are legally obliged to return old batteries and accumulators after use. You can do this free of charge at a collection point in your area. Addresses of suitable collection points can be obtained from your city or local government.
3. Batteries may contain substances that are harmful to the environment and human health. The separate collection and recycling of old batteries and accumulators is intended to avoid negative effects on the environment and human health.
4. Please avoid the generation of waste from old batteries as far as possible, e.g. Please avoid littering public spaces by not carelessly leaving batteries or electrical and electronic equipment containing batteries.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/30/EU, Directive 2011/65/EU, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address: <https://renpho.eu/pages/compliance-data>