Dreame Chrona Pro Photon Anti-Aging Skincare Mask

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Warning

To prevent electric shock, fire, or other accidents, please read this user manual carefully before use and retain it for future reference.

- Before using the LED light therapy mask, you must perform a sensitivity test (please refer to the Photosensitivity Test section). Do not use this product if you experience an allergic reaction when using a LED light therapy mask. This device is intended for use in a home healthcare environment and may be operated by the user.
- Do not use the LED light therapy mask in any manner not explicitly specified in this manual. Any use of the device other than what is explicitly specified in this manual has not been tested, and the associated risks are unknown. This device is not intended for use by individuals with reduced physical, sensory, or mental capabilities, or a lack of experience and knowledge, unless they are supervised or instructed by a responsible person to ensure safe use and to avoid risks of fire or burns.
- If you feel the device is overheating or your skin feels uncomfortable due to the heat, immediately turn off and remove the LED light therapy mask.
- This product operates wirelessly during use. Only connect the charging cable when charging.
- During charging, the charging cable may pose a risk of entanglement or strangulation. Always keep the device and charging cable out of reach of children.
- Remove makeup before using the LED light therapy mask. Certain cosmetic ingredients may cause photosensitivity or reduce the effectiveness of the treatment. Using the mask without thoroughly removing makeup may cause skin reactions or reduce the effectiveness of the treatment. This device is not intended for use by children. Ensure close supervision to prevent children from playing with the device to avoid the risk of fire and burns.
- Do not use the LED light therapy mask near water sources or in environments where it may be splashed with water. Do not use the mask if you suspect it has come into contact with water. Avoid strong impacts to prevent damage to the light modules.

- Keep the LED light therapy mask away from hot surfaces. If necessary, contact the manufacturer for assistance with installation, use, or maintenance, or to report abnormal operations or incidents.
- Do not use the LED light therapy mask if the mask, controller, or USB Type-C charging cable is damaged. Do not leave the device unattended when powered on to avoid fire, burns, or other hazards.
- When the product is taken from minimum or maximum storage temperature conditions, allow it to rest for 2 hours before use.
- Do not stare directly into the light during use or allow it to shine directly into the eyes.
- Do not use the LED light therapy mask near devices using wireless technology.

Precautions

- This device contains silicone, plastic components, and hook-and-loop fastener materials. If you have a known allergy to any of these materials, discontinue use and consult a healthcare professional.
- This device features a memory function. Each time it is turned on, it automatically recalls and restores the last used treatment mode, intensity level, and treatment duration.
- LED light therapy mask cannot be used while charging.
- Before each use, clean the mask with the device powered off.
- Do not use any cleaning products other than those mentioned in the cleaning instructions. Using other cleaning products may damage the LED light therapy mask.
- Store the LED light therapy mask in a dry environment at room temperature, away from direct sunlight.
- When the product is taken from minimum or maximum storage temperature conditions, allow it to rest for 2 hours before use.
- Do not stare directly into the light during use or allow it to shine directly into the eyes.


Electrical Safety

- Do not damage, modify, bend, twist, or heat the USB Type-C charging cable. A damaged charging cable may cause fire, electric shock, or malfunction.
- Do not place or operate the device near water or other liquids. Otherwise, it may cause fire or electric shock.
- This product uses a magnetic detachable battery pack. Please fully charge the battery pack before using the device for the first time.
- Only use the dedicated battery pack provided with the device.
- Do not attempt to open, modify, replace, or use unapproved battery packs. Improper handling may cause overheating, leakage, or explosion.
- Do not dispose of batteries improperly. Please follow local regulations for proper disposal.
- Only use a DC 5V $\equiv$ 1A charger to charge this device (maximum output power 3.5W).
- Do not use chargers of unknown origin, inferior quality, or with an output power exceeding 30W, otherwise it may cause overheating or damage to the product. To ensure maximum safety, it is recommended to use a charger that complies with IEC 60601-1:2005 + AMD2:2020 standard, and a charger with the following specifications: Input: 100-240VAC, Output: DC 5V $\equiv$ 1.0A.

The product body is marked with "Do not use while bathing or showering".

EN Product Overview

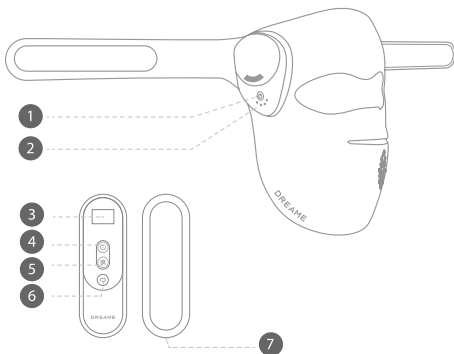1 Power Button

2 Battery Indicator

4 green lights on: 100% battery (fully charged)
3 green lights on: Approx. 75% battery remaining
2 green lights on: Approx. 50% battery remaining
1 green light on: Approx. 25% battery remaining
Red light flashing: Battery below 10% (low battery)%

3 Display

Used to view the selected mode/intensity/timer

4 Timer Button

5 minutes/10 minutes/15 minutes

5 Intensity Level Button

Level I (low intensity)
Level II (medium intensity)
Level III (high intensity)

6 Power/Mode Button

Five modes: Firming Mode, Brighten Mode, Spot Fading Mode, Soothing Mode, Multi-effect Mode

7 Type-C Port

Note: All product images, accessories, and user interface illustrations in this manual are for reference only. Due to product updates and upgrades, the actual product may vary slightly from the illustrations. Please refer to the actual product.

How to Use

Setting Up the Device

1. Take out the device and all included accessories.
2. This product is equipped with two magnetic detachable batteries (1000mAh each, 2000mAh in total). Before first use, fully charge both batteries and the controller. For specific operations, please refer to the "Charging Instructions" section.

Photosensitivity Test

Please follow the steps below to test if you are sensitive to the treatment light:


1. Press the power button to turn on the device. Use the controller to turn the device on/off and select the desired treatment mode.
2. Initial test: Place the beauty mask closely against a small area of skin on your forearm and complete a full skincare session (3 minutes) in the selected mode.
3. Switch to another mode and perform a second session (3 minutes) on the same skin area.
4. After completing the test, wait 24 hours. During this time, the tested area should not show any adverse reactions.

Adverse reactions include pain or discomfort at the treatment site, persistent redness, and other abnormal skin reactions.

EN Charging Instructions

This device is equipped with two detachable 1000mAh battery packs and a wireless controller (with a built-in 280mAh battery).

Before first use, please refer to the illustration below:

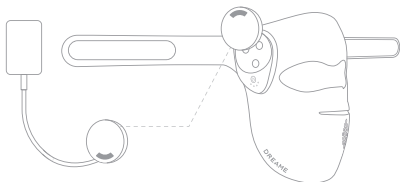Battery Pack Charging Instructions

- Remove the battery pack from the mask before charging.
- Each battery pack takes approximately 2 hours to fully charge.
- Actual charging time may vary depending on the remaining battery level.

Controller Charging Instructions

- This wireless controller features a built-in 280mAh rechargeable battery.
- It takes approximately 2 hours to fully charge.
- Actual charging time may vary depending on the remaining battery level.

Lithium Battery Maintenance Tips

- Avoid completely draining the battery before recharging.
- Avoid overcharging the battery; do not charge for more than 8 hours.
- If the device is not used for an extended period, it is recommended to charge the battery at least once every three months.

Getting Started

Before Using the Device

1. Cleanse your face (remove makeup first if applicable).
2. This product can be used morning and evening, immediately after cleansing or incorporated into your daily skincare routine.
3. Ensure the device is fully charged before use.

Procedure

1. Short press the power button to turn on the mask (it can be turned on when the green light flashes).
2. Adjust the headband appropriately and put on the mask, ensuring the inner side fits snugly against your face.

Headband Adjustment Instruction

1. When using the device in a sitting position, adjust the headband to minimize pressure on the eye sockets while ensuring the silicone eye pads fit snugly around the eyes to secure the product.
2. If the silicone eye pads put too much pressure on your eye sockets, slightly loosen the headband on both sides.
3. If light leaks through a gap between the eye pads and your eye sockets, tighten the headband on both sides.

4. If the mask slips easily or fails to stay in place, slightly tighten the headband.
5. You can select the desired treatment mode using the controller.
6. If you experience any discomfort while wearing the mask, try adjusting the headband or removing the mask to rest until you feel comfortable.
7. The device will automatically shut off after the set mode time ends.

Precautions

- The controller must be used within 2 meters of the mask.
- Avoid strenuous activities while wearing the mask.
- During use, you can lie down or lean back on a sofa or bed to relax.

Controller Instructions

1. Connect the included charging cable to the Type-C port on the controller.
2. Plug the USB end of the charging cable into a 5V 1A power adapter for charging.
3. Charging status indicators:
 - "88" flashing → Charging in progress
 - "88" solid → Charging complete(Unplug the charging cable before use)
4. Low battery alert: When the controller battery is low, "88" will flash continuously until the screen turns off. Please charge it promptly.

Since the LED mask uses a wireless RF controller (the controller emits RF signals, and the main unit receives them), users need to pair the controller with the LED mask before first use.

The pairing steps are as follows: Press and hold the power button on the LED mask for 5 seconds to enter pairing mode. The indicator light on the mask will switch from green to red, and then flash rapidly between green and red.

Power On/Off Operation

Mode Selection

Press the "Power/Mode Button" (which serves as both the power switch and mode selector) to switch between the following treatment modes:

Mode 1: Firming Mode Red Light (630nm) + Near-Infrared Light (850nm) + Near-Infrared Light (1070nm)

Benefits: Targeting early signs of aging and loss of firmness, the triple light energy works synergistically to penetrate deep into the skin and stimulate collagen renewal. It significantly improves elasticity and firmness, accelerates skin rejuvenation, and restores a smooth, plump, and youthful radiance.

The 1070nm wavelength offers deep penetration, reaching the deep dermis, subcutaneous tissue, and even the superficial muscle layer. Compared with shallower light waves, it targets deeper skin structures to help improve fine lines and sagging while enhancing skin firmness and elasticity.

Mode 2: Brighten Mode Yellow Light (590nm)

Benefits: Soothes and revitalizes the skin, making it ideal for skin recovery after sun exposure or fatigue. Helps brighten the complexion, even out skin texture, and restore a healthy glow.

Mode 3: Spot Fading Mode Green Light (530nm)

Benefits: Calms and balances sensitive or stressed skin, reduces redness, and promotes an even skin tone. Leaves the skin naturally balanced and refreshed, especially suitable for reactive or sensitive skin.

Mode 4: Soothing Mode Blue Light (465nm)

Benefits: Commonly used in skincare routines, it helps soothe visible redness, reduce the appearance of shine, and support a clearer-looking complexion. It keeps the skin refreshed and balanced, helping to promote a healthier, more refined skin texture.

Mode 5: Multi-effect Mode (Yellow Light 590nm + Green Light 530nm + Red Light 630nm + Near-Infrared Light 850nm + Near-Infrared Light 1070nm)

Wavelength Combination: Yellow Light (590nm) + Green Light (530nm) + Red Light (630nm) + Near-Infrared Light (850nm + 1070nm)

Benefits: Integrates multi-wavelength light energy to synergistically achieve four major skincare goals: wrinkle reduction, brightening, dark spot reduction, and purifying. Comprehensively revitalizes the skin from within, evens skin tone, reduces hyperpigmentation, and improves overall skin texture. As an all-in-one daily facial treatment, it provides both repairing and preventive benefits for a lasting youthful glow.

Cold Compress Instructions

You can choose an appropriate cold compress method based on your personal needs and skin condition:

I. Prepare the Cold Pack

Gentle Cold Compress

For soothing relaxation and fatigue relief, lay the cold pack flat in the refrigerator (0°C -10°C) for about 15 minutes before use. The temperature of the cold pack will be approximately 20°C .

* Note: For daily use, it is recommended to keep the cold pack in the storage bag and store it in the refrigerator compartment.

Intensive Cold Compress

For rapid cooling or an enhanced cooling effect, lay the cold pack flat in the freezer compartment (-15°C to -20°C) for about 10 minutes before use. The temperature of the cold Pack will be approximately 15°C.

Storage Tips

It is recommended to keep the cold pack in the refrigerator compartment for daily use and only place it in the freezer compartment when intensive cooling is required. Do not leave it in the freezer for more than 15 minutes, otherwise ice crystals may form.

II. How to Use

Remove the cold pack from the refrigerator, attach its magnetic side to the magnetic area of the mask, and then apply it to the desired area.

The recommended duration for each cold compress session is 15-20 minutes, which can be adjusted based on your personal tolerance. If continuous cold therapy is needed, please wait at least 1-2 hours before the next session.

III. Functional Description of the Two Cold Packs

This product includes two accessories: a full-face cold pack and a targeted cold pack, to meet the cold therapy needs of different areas.

Full-face cold pack


Suitable for cooling large areas of the face, helping to alleviate overall facial heat, redness, or slight swelling. It is ideal for uniform cooling of the entire face.

Note: The full-face cold pack comes with a strap. Due to individual differences in facial contours, if it does not fit perfectly against your face during application, you can use the strap to secure it.


Targeted cold pack (for areas under the eye)

Specially designed for the eye area, it fits the under-eye contours to relieve puffiness, dark circles, or fatigue. It is particularly suitable for soothing cold therapy on the eyes after staying up late, eye strain, or cosmetic treatments.


Cold Compress Instructions

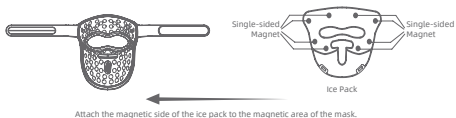Full-face cold pack

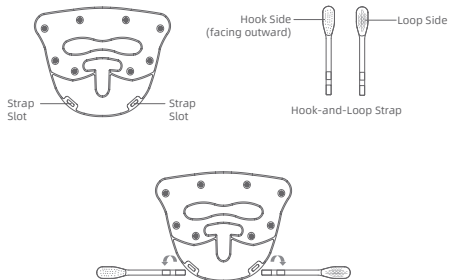Targeted cold pack

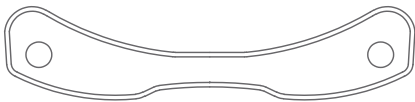Remove the cold pack from the refrigerator, attach its magnetic side to the magnetic area of the mask, and then apply it to the desired area. The recommended duration for each cold compress session is **15-20 minutes**, which can be adjusted based on your personal tolerance. If continuous cold therapy is needed, please wait at least 1-2 hours before the next session.

How to Use

Use with the beauty mask

After using the LED beauty mask, you can choose the full-face or targeted cold pack for a cold compress session based on your skin condition to help soothe and calm the skin.

Standalone use

The cold pack can also be used independently without the mask. It is suitable for daily scenarios such as relieving eye fatigue, reducing facial puffiness, or cooling slight bumps, offering flexibility and convenience.

* You may also detach the cold pack to use the mask alone for light therapy.

Cold Pack Precautions

Discontinue use immediately if you experience numbness, tingling, or if your skin turns white while using the cold pack.

If you have broken skin, eczema, or blood circulation disorders, please use with caution or under the guidance of a doctor.

Children, the elderly, or individuals with reduced sensation should be accompanied by an adult and closely monitored when using the cold pack.

EN Changing the Intensity Level

The mask offers the following intensity levels:

Intensity 1 (low)

Intensity 2 (medium)

Intensity 3 (high)

Timer Settings

Press the timer button to select the following treatment durations:

· 5 minutes

· 10 minutes

· 15 minutes

Memory Function Description

For the first use, the device is preset to Mode 1 (revitalizing red light), Intensity 3 (high), and a 10-minute timer. After each use, the device automatically saves the previously used treatment mode, intensity level, and timer setting.

Notes

1. The mask will automatically turn off after the set treatment time ends.
2. When the controller is turned off, the mask will also turn off.
3. If you select a treatment mode via the controller and later press the mask buttons again, the mask will default to a 10-minute timer. In this case, the timer displays on the mask and the controller will not match. Therefore, we recommend starting the light therapy immediately after selecting a mode with the controller.

Cleaning and Maintenance

1. Please store the mask in a dry place at room temperature. We recommend using the provided storage bag.
2. Do not store the product in a dusty or humid environment, or under direct sunlight.
3. Do not disassemble or modify the internal components.
4. Before each use, turn off the device and wipe the contact surface of the mask according to the cleaning instructions below.
5. Avoid excessive bending or twisting of the silicone part of the mask.
6. Do not apply cosmetics directly onto the mask surface.
7. Avoid contact with oils, lotions, and perfumes during use.

Cleaning Instructions

1. The mask body is made of silicone. After each use, gently wipe the surface with a soft, lint-free cloth dampened with water.
2. For deep cleaning, gently wipe the mask with a mild, alcohol-free, and fragrance-free neutral soap solution diluted with water, then wipe it clean with a clean, damp cloth.

Daily Cleaning of the Cold Pack

1. After each use, gently wipe the surface of the pack with a soft, damp cloth or sponge moistened with a small amount of neutral detergent (such as dish soap) to remove sweat or dirt.
2. Thoroughly wipe off any detergent residue with clean water to prevent chemicals from corroding the PVC material.
3. Hang the cold pack or lay it flat in a cool, well-ventilated area to air dry. Do not expose it to direct sunlight, tumble dry, or blow dry with hot air, as this may accelerate aging.

Cold Pack Storage and Pre-Cooling

Recommendations

Daily storage: After cleaning and air drying, store the cold pack in the refrigerator compartment (not the freezer). **Do not leave it at room temperature outside the refrigerator,** as this can cause the internal cooling gel to evaporate quickly and draw in air, creating bubbles in the cold pack and affecting its performance.

Precautions

Avoid contact with sharp objects to prevent punctures and leaks.

Do not wash with boiling water or heat in a microwave, as the PVC material cannot withstand high temperatures.

Regularly check the seal of the pack. If any damage is found, stop using immediately and replace it.

Disposal Instructions

Do not dispose of any part of this product as regular household waste.

This product and its controller contain lithium-ion polymer batteries. To protect the environment, at the end of the product's life, please:

Return the fully packaged device to the manufacturer or its authorized representative, or strictly follow your local electronic waste disposal regulations for safe recycling.

FAQs and Troubleshooting

EN

Issue	Solution
The controller screen is off and it is unresponsive and cannot control the mask.	- Check if the controller's battery is depleted. Connect the charging cable to the Type-C port to charge it. - If it still does not work after charging, please contact the manufacturer or customer service center for inspection.
The controller screen is off, but it can control the mask.	Follow the pairing steps to reconnect the controller to the mask. - Press and hold the power button on the mask for 5 seconds. The indicator light will first turn green, and then flash red, indicating that the device is in pairing mode. - Press the power button on the controller to complete pairing.
The mask cannot be charged.	- Try a different charging cable. If the charging cable is causing this issue, contact us for a replacement. If the issue persists, please contact us to replace the device. - Check the power adapter and voltage. Please check if the power adapter connected to the device is receiving power, and ensure that the voltage indicated on the device matches that of your local power supply.
Bubbles are found in the cold pack or the cold pack weight is reduced by 50%.	- Replace it with a new cold pack.

Product Specifications

Product Name	Dreame Chrona Pro Photon Anti-Aging Skincare Mask
Model Name	AB01P
Battery Specifications	Type: 3.7V / 1000mAh lithium polymer battery (each) Quantity: 2 detachable battery packs (total capacity: 2000mAh)
Product Size	Mask: 634*233.8*5.5mm Controller: 133.6*42.3*19mm
Product Net Weight	Mask: About 338g Controller: About 72g Full-face cold pack: About 155g Targeted cold pack: About 18g
Charging Time	About 2 hours
Battery Usage Time	One fully charged battery: About 60 minutes Two batteries combined: About 120 minutes in total

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Manufacturer: Dreame Trade (Tianjin) Co., Ltd.